

Multiple Applications

- Pharmaceutical Research and Development
- Drug Stability and Shelf Life Testing
- Plasma and Blood Storage
- Clinical Research
- Environment Simulation Testing
- Working Laboratories
- Biological Research
- Electronic Stress Testing


Mini-Room Walk-in available as...

- Cold room, freezer, combination cold room/freezer or warm room models
- Indoor or outdoor models
- Sizes from 4' x 5' to 8' x 14'
- Three Heights: 6'7" & 7'7" with floor, 7'4" floorless
- Capacity from 77 cubic feet to 1,000 cubic feet
- Wall or ceiling mounted Capsule Pak™ Refrigeration Systems
- Temperatures: +4°C, -20°C, -30°C, +2°C to +55°C



The Mini-Room Walk-in offers you...

- Installation that is fast and easy
- A full 4" panel foamed-in-place with CFC and HCFC free HFC-245fa polyurethane insulation
- Deadbolt locking handle for security
- Nor-Lake Scientific's exclusive Capsule Pak™ refrigeration system that features a flush coil. 100% usable interior space
- Automatic condensate vaporizer. No plumbing required on indoor units
- Expanded Mini-Room Cold Rooms
 - 6' x 14'
 - 8' x 12'
 - 8' x 14'
 - 8' x 20' combination
- Smaller Mini-Room Walk-ins
 - 3'6" x 6'0" x 6'7" and
 - 3'6" x 7'0" x 6'7"
- Similar floor space as 3 door laboratory reach-in with 20%-30% more capacity

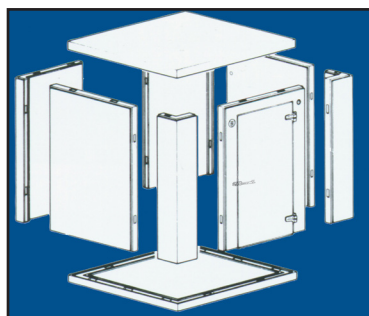

Modular Assembly

Capsule Pak™ Refrigeration Systems

- Capsule Pak™ Refrigeration Systems consist of a unitized system which is factory assembled, wired, charged, tested and fully equipped for insertion into a factory prepared walk-in opening. The Capsule Pak™ System has a flush coil which keeps all components outside, allowing more storage inside.
- Wall or ceiling mounted systems
- Indoor or outdoor models
- Temperatures: +4°C, -20°C, -30°C, +2°C to +55°C
- Optional Remote Capsule Pak™ Systems
- UL & C-UL listed

Shelving Kits

Specialty designed Shelving Kits are available for all Mini-Room walk-ins.

- Nor-Lake Scientific Chrome Kote™ Shelving Systems
- Chrome Kote™ shelving available as standard when included
- Quick, easy assembly
- Rigid construction
- Stronger shelves
- Flexible shelf adjustment- 1" increments
- Easy to use
- Nor-Lake shelving utilizes a support clip to eliminate center posts
- The Nor-Lake Klip™ provides maximum space utilization for storage leaving the door opening clear of end standards
- Special epoxy coated finish Green Kote™ shelving provides extra protection from acidic products.

Contact your sales representative for additional details.

NORLAKE SCIENTIFIC

727 Second Street
Hudson, Wisconsin 54016
800-477-5253
715-386-2323
715-386-4290 FAX
www.norlakescientific.com
© 2014 NOR-LAKE, INC.
Revision Date: 06/14
Part Number: 089730

ISO 9001:2008
Registered to ISO 9001:2008
File No. 10001816

Standex
Food Service Equipment Group

